



Molift Mover 180

A lightweight and compact hoist, suitable in home care and institutional environments

Easy to manoeuvre

The Molift Mover 180 features a lightweight aluminium construction, making it easy to manoeuvre. Its large hoisting interval of 27–163 cm allows for hoisting from both high surfaces and the floor. Both 2-point and 4-point sling bars are available and can easily be changed without tools thanks to the quick release coupling,

Equipped with 75 mm castors, it provides excellent manoeuvrability in tight spaces. For movement over uneven surfaces, 100 mm castors are available, enhancing stability and ease of use.



molift[®]



Hand control

Robust and washable, Indicators for battery, overload, and service.



Unique push bar

The push bar ensures ergonomic, natural grip with multiple options.



Sling bar options

2-point and a 4-point sling bars are available in several sizes.



Easy to manoeuvre

Lightweight and compact design makes it easy to manoeuvre.

Features

Ergonomic handling

The Molift Mover 180 is equipped with Molift's unique push bar that provides multiple and flexible grip possibilities for improved caregiver ergonomics.

Life expectancy

The hoist has an expected lifetime of 30,000 cycles or 15 years, whichever comes first.

Battery

The battery features advanced Li-Ion technology and a safe, sealed design makes it washable (IPX6). It weighs only 0.9 kg, and can be replaced with a single hand.

The battery can be charged directly on the hoist, and an external charger is also available. The charging time is only 3.5 hours. Fully charged it delivers 160 hoists at 75 kg per, and is compatible with all Molift mobile hoists.

With a life expectancy of 1,700 charging cycles and a 5-year warranty, the battery is built to last.



Speed boost

The speed boost feature enables increased speed, when the actuator is unloaded.

Hand control


The hand control is robust and well-proven. It features a battery level indicator, overload indicator, and service indicator to ensure optimal performance and safety.

Designed for durability, the cable can be easily replaced if damaged, and the hand control is washable (IPX6).

Molift Service Tool

Molift Mover 180, like all Molift hoists, is equipped with a unique service notification system that notifies when service is required. The integrated service software can also display data about hoist usage and can assist service technicians in troubleshooting.

Additionally, the hoist is equipped with Bluetooth connectivity, enabling wireless connection to Molift Service Tool.

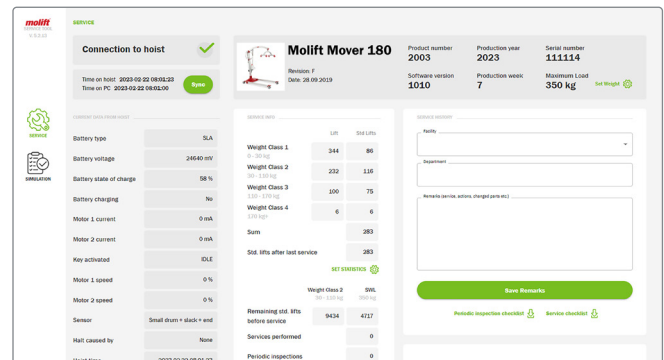

Product Molift Mover 180 Molift Mover 180 - excluding sling bar (ordered separately)		Item no. M1801 M1802	Dimensions Outside measurement with closed legs 670 mm Outside measurement with open legs 1145 mm Inside measurement with closed legs 555 mm Inside measurement with open legs 1025 mm
Material Aluminium, steel, plastic			Technical Battery Li-Ion 25.2V, 2.85Ah Battery charge time 3.5 hours Hoisting speed 46 mm/sec Rated performance 160 lifts/75 kg
Weight Safe working load (SWL) 180 kg Weight of unit (excl. battery) 31 kg			Cleaning Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution
Dimensions Height 1940 mm Length 1260 mm Leg height 105 mm Castor size 75/100 mm Hoisting interval 1365 mm Hoisting min-max with 2-point sling bar 360 - 1730 mm Hoisting min-max with 4-point sling bar 265 - 1630 mm			

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com