



Molift Mover 300

When a high safe working load is needed

With its high hoisting capacity of 300 kg, Molift Mover 300 is the ideal solution for hospitals and care facilities with bariatric clients who have frequent and everyday moving and handling needs. With a total weight of only 36 kg the Molift Mover 300 is the lightest in its category. The unique construction and wheel positioning provides excellent manoeuvrability.

Molift Mover 300 hoists easily from the floor, reaching down to 30 cm at the lowest. Excellent maximum hoisting height at 171 cm provides good free hoisting space and possibility for gait training. The hoists battery can be replaced with a single hand and is washable (IPX6).



molift[®]



Hand control

Robust and washable. Indicators for battery, overload, and service.



Unique push bar

The push bar ensures ergonomic, natural grip with multiple options.



Quick-release function

The quick-release function makes it easy to change to the sling bar.



Sling bar options

2-point and a 4-point sling bars are available in several sizes.

Features

Ergonomic handling

The Molift Mover 300 is equipped with Molift's unique push bar that provides multiple and flexible grip possibilities for improved caregiver ergonomics.

Life expectancy

The hoist has an expected lifetime of 30,000 cycles or 15 years, whichever comes first.

Battery

The battery features advanced Li-Ion technology and a safe, sealed design makes it washable (IPX6). It weighs only 0.9 kg, and can be replaced with a single hand.

The battery can be charged directly on the hoist, and an external charger is also available. The charging time is only 3.5 hours. Fully charged it delivers 110 hoists at 75 kg, per full charge, and is compatible with all Molift Mover hoists.

With a life expectancy of 1.700 charging cycles and a 5-year warranty, the battery is built to last.

Speed boost

The hoist comes with a speed boost feature that enables increased speed, when the actuator is unloaded.


Hand Control

The hand control is robust and well-proven. It features indicators for battery level, overload, and service, to ensure optimal performance and safety. Designed for durability, the cable can be easily replaced if damaged, and the hand control is washable (IPX6).

Molift Service Tool

Molift Mover 300, like all Molift hoists, is equipped with a unique service notification system that notifies when service is required. The integrated service software can also retrieve data about hoist usage and can assist service technicians in troubleshooting.

Additionally, the hoist has Bluetooth connectivity, enabling wireless connection to Molift Service Tool.

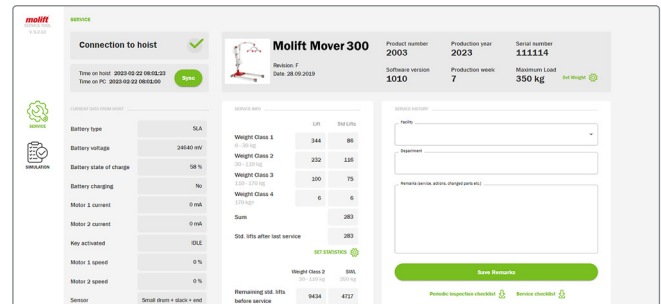

Molift Mover 300		Item no.	Dimensions
Molift Mover 300		M1501	Castor size 100 mm
- excluding sling bar			Hoisting interval 1390 mm
Molift Mover 300		M1502	Hoisting min-max with 2-point sling bar 400 - 1815 mm
- excluding sling bar, low base			Hoisting min-max with 4-point sling bar 300 - 1715 mm
Molift Mover 300		M1503	Outside measurement with closed legs 690 mm
- excluding sling bar, incl. support arms			Outside measurement with open legs 1205 mm
Molift Mover 300		M1504	Inside measurement with closed legs 550 mm
- excluding sling bar, incl. support arms, low base			Inside measurement with open legs 1065 mm
Material			Technical
Aluminium, steel, plastic			Battery Li-Ion 25.2V, 2.85Ah
Weight			Battery charge time 3.5 hours
Safe working load (SWL)		300 kg	Hoisting speed 32.6 mm/sec
Weight of unit (excl. battery) M1501 & M1502		36 kg	Rated performance 110 lifts/75 kg/
Weight of unit (excl. battery) M1503 & M1504		42.3 kg	
Dimensions			Cleaning
Height		2025 mm	
Length		1285 mm	
Leg height		115 mm	
Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution			

For the intended purpose of the products and latest information, please visit www.etac.com

Etac

+46 371 587 330

info@etac.se

www.etac.com